



Iberose[®] - A Leading Brand Of Electrophoresis Grade Agarose

Iberose is the brand name used by Web Scientific for a range of related products which cover a variety of applications from routine genomic analysis: analysis of DNA and RNA fragments, cloning, PCR, blotting, etc (General Purpose Iberose) to isolation of DNA fragments from gels (LMP Iberose) to molecular screening applications in conjunction with high throughput electrophoresis systems (MS8 Iberose). All our products are manufactured to the highest standards by one of the world's leading producers of agarose, so quality together with research grade specifications are things you can be sure of.

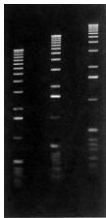


Iberose[®]

Iberose is a general purpose agarose suitable for routine analysis of nucleic acids. Ideal for routine rapid separation of DNA and RNA fragments as well as PCR products, the preparation of plasmids, and for screening, cloning and blotting techniques.

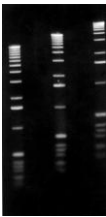

Iberose[®] LMP

Low Melting Point (LMP) Iberose allows for the recovery of undamaged nucleic acids below the denaturation temperature. Recommended for $\geq 1,000$ bp fragment separation in analytical electrophoresis. The low melting temperature of this agarose (65°C) allows the gel to be melted without damaging the DNA, so that fragments can be recovered. In-gel Manipulations can also be performed without prior extraction of the DNA from the gel. DNase/RNase activity: none detected.

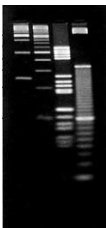

Iberose[®] MS8

Molecular screening Iberose (MS8) can be used to improve resolution of small DNA fragments and PCR products. Recommended for analytical gels for DNA $\leq 1,200$ bp and especially for PCR products. DNase/RNase activity and DNA binding: none detected.

Ordering Information

Prod. Code	Description	Pack Size	Price
AGR-100	Iberose - Agarose for electrophoresis	100g	£34
AGR-500	Iberose - Agarose for electrophoresis	500g	£150
AGR-2000	Iberose - Agarose for electrophoresis	4 x 500g	£462
LMP-100	Iberose LMP - Low Melting Point Agarose	100g	£198
MS8-100	Iberose MS8 - Molecular Screening Agarose	100g	£165

***Note: Other pack sizes of Iberose are available. Please enquire for further details**